

U-10 **DILUTE** INSULIN (10 UNITS/ML) CALCULATOR HELP

The CHLA **DILUTE** Insulin Calculator will accurately determine the correct dose of **dilute (10 units/mL)** rapid-acting insulin to be administered based on a specific patient order, the measured blood glucose, and (IF carb-counting) the grams of carbohydrates eaten.

Dilute Insulin is generally used in very small patients who are very sensitive to insulin. Because Dilute insulin is only $1/10^{\text{th}}$ the concentration of standard U-100 insulin, **doses are expressed in Syringe Units** based on markings on an insulin syringe.

Insulin/Carbohydrate Ratio: (ABSENT IF CORRECTION-ONLY ORDER)

- 1 Syringe Unit of insulin covers **X** grams of carbohydrates eaten
- Allowed range: 1.5 – 40 gm: 1 Syringe Unit
 - ✓ ordered value; click to confirm → ✓
 - ✗ indicates an ordered value that should be clarified

Correction Factor (CF) -- also called Insulin Sensitivity Factor (ISF):

- 1 Syringe Unit of insulin reduces a Blood Glucose elevated over Target Blood Glucose by **X** mg/dL
- Allowed range: 10 – 200 mg/dL: 1 Syringe Unit
 - ✓ ordered value; click to confirm → ✓
 - ✗ indicates an ordered value should be clarified

Target Blood Glucose:

- blood glucose target above which correction must be applied
- Allowed range: 130 – 250 mg/dL
 - “Usual” Target BG ($\sim 10^{\text{th}}$ – 90^{th} %ile): 150 – 200
 - ✓ ordered value; click to confirm → ✓
 - ⚠ indicates an extremely low or high value that should be clarified
 - LOW: < 130 | HIGH: > 250

***The calculator is not enabled for input of carbs eaten or measured blood glucose until all above parameters are confirmed by clicking the ✓ or ⚠**

Rapid-acting insulin in past 2 hrs?:

- *If yes, will not include any correction dose for high blood glucose in total dose calculation in order to prevent insulin stacking for insulin on board*
- Calculated dose will include carb coverage only
 - Defaults to YES if MAR shows rapid-acting insulin charted within past 2 hrs
 - Defaults to NO if recent dose not detected
 - No recent dose
 - Dose given but not charted (home, pump, other)
- response may be changed based on circumstance


Carbohydrates Eaten: (ABSENT IF CORRECTION-ONLY ORDER)


- *Enter the grams of carbohydrates eaten*
- calculator will accept any value from 0 – 99 gm
-  unacceptable entry which must be corrected

Blood Glucose for Calculation:

- *Enter the current glucose measurement*
- calculator will accept any value from 1 – 2,000 mg/dL
-  BG > 600 mg/dL; verify and correct if data entry was incorrect
-  BG < 70 mg/dL: hypoglycemia
 - STOP – Calculator will not calculate dose
 - DO NOT GIVE INSULIN
 - requires rescue intervention for low BG
-  unacceptable entry which must be corrected

To complete calculation of dose to administer:

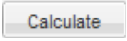- *Assure all parameters and patient values are valid*
- Click  button
- Dose will be displayed, **rounded** to closest Syringe Units, here:

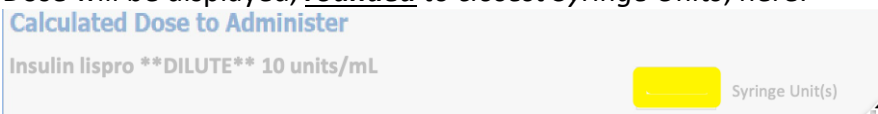

Calculated Dose to Administer

Insulin lispro **DILUTE** 10 units/mL

 Syringe Unit(s)

- Transfer calculated dose to MAR after independent 2nd witness check and dose has been administered
- **REMEMBER WHEN CHARTING**

- *Change time from default to time actually administered*
- *Change dose from “**1 dose**” to “**X Syringe Units**” actually administered*